
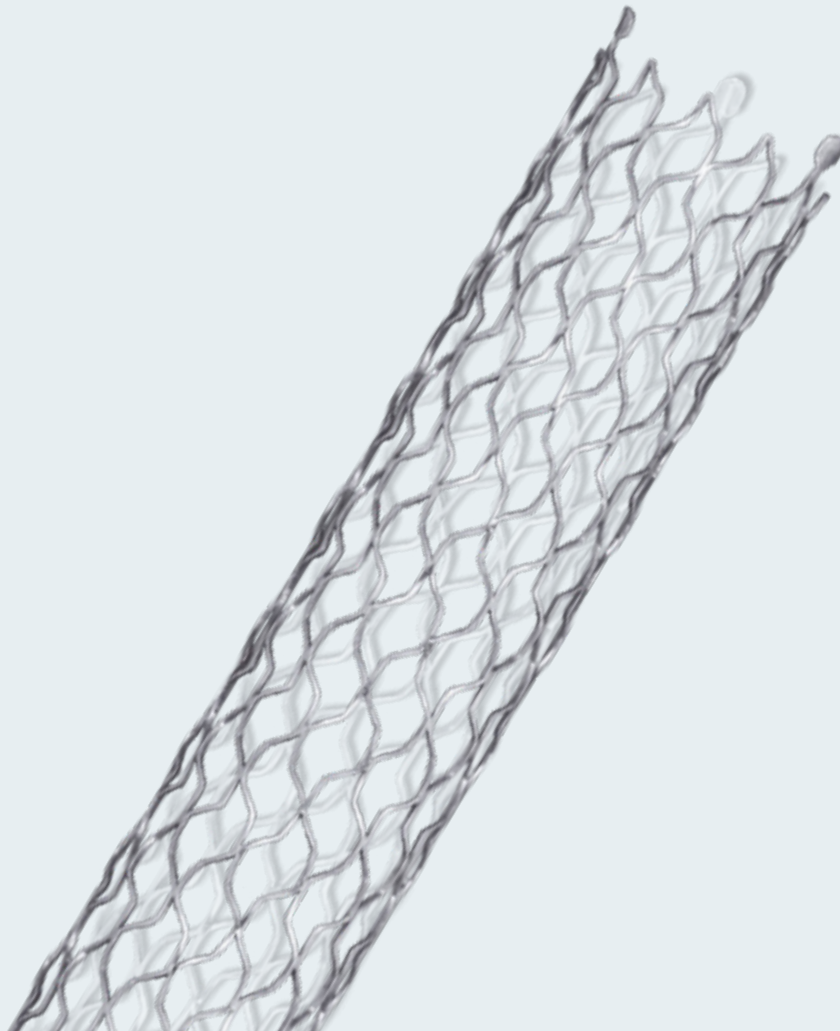
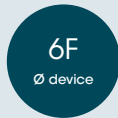


sinus-Repo- Visual 6F

REPOSITIONABLE & POWERFUL -

**A DEDICATED STENT DESIGN FOR STRAIGHT
VESSEL SEGMENTS IN THE ILIAC ARTERIES.**


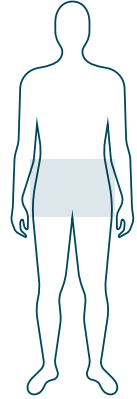


sinus-Repo-Visual 6F

REPOSITIONABLE AND POWERFUL

INTENDED USE

The sinus-Repo-Visual 6F vascular stent system is intended for the endovascular treatment of atherosclerotic lesions or obstructions in the iliac arteries.

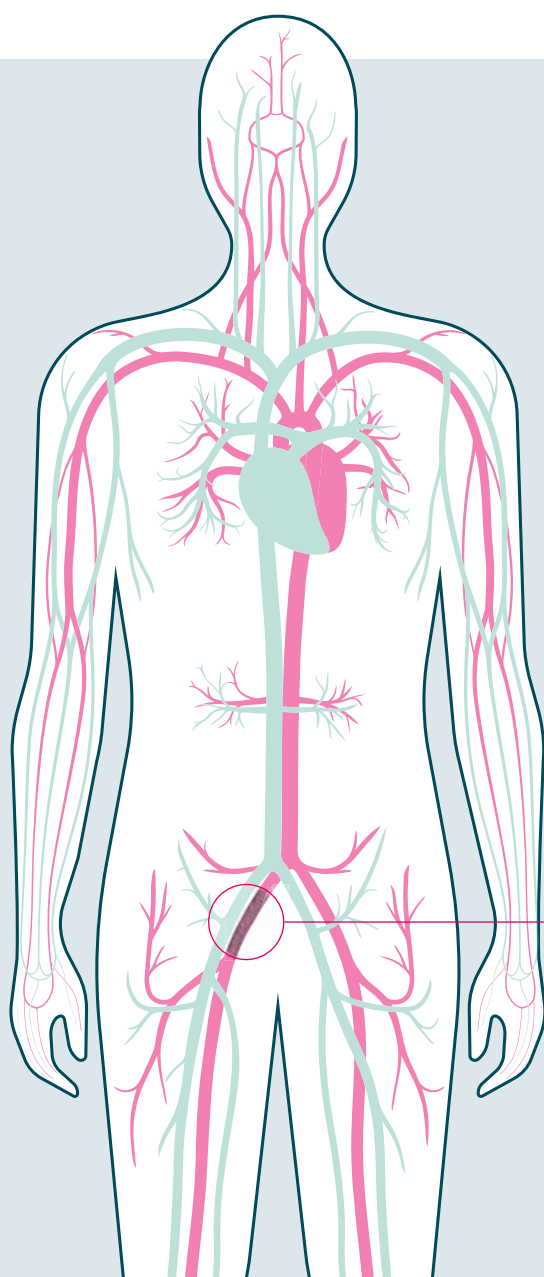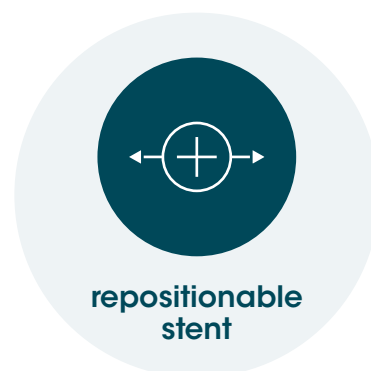

closed-cell
design




radiopaque
markers



anti-jump
technique™



electro-
polishing

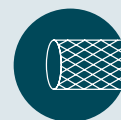

braided
sheath




atraumatic
soft-tip



REPOSITIONABLE DESIGN

The sinus-Repo-Visual-6F vascular stent system has been specifically designed to enable a repositioning of the stent in the target lesion - this special feature facilitates a most precise stent placement.

The repositioning mechanism is supported by 2 markers on the delivery system that indicate the area in which repositioning is enabled. As long as the red marker (point of no return) is visible, repositioning of the stent is possible.

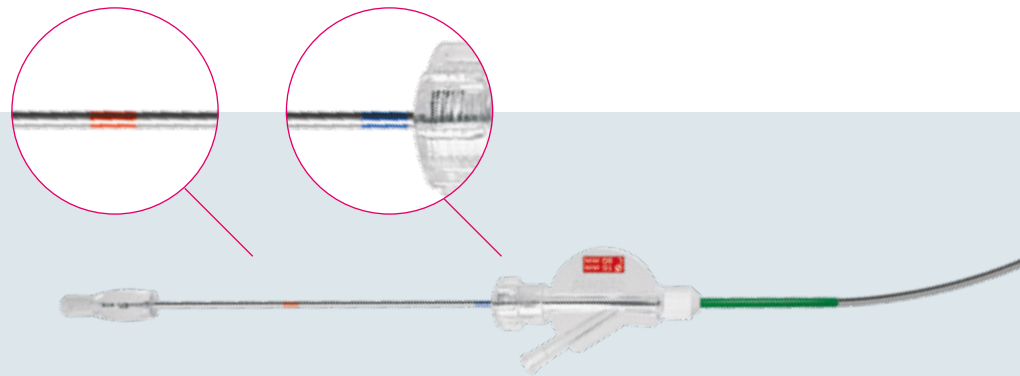

FEATURES


6F delivery system adapted to 0.035 inch guide wires

Repositionability feature for an even more precise stent placement increasing the chances of successful procedure

Maximum radial force for straight vessel segments and highly calcified stenosis

Closed-cell stent design for a high lesion coverage

LENGTHS

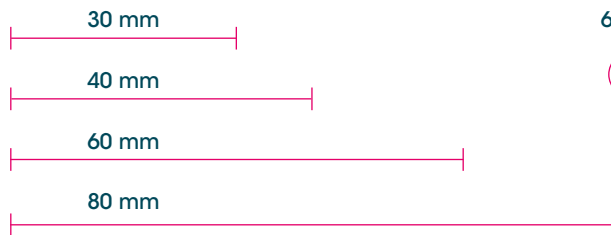

DIAMETERS



ORDER CODES - 75CM APPLICATION DEVICE

Stent Ø (mm)	Stent length (mm)			
	30	40	60	80
6	6306-6030	6306-6040	6306-6060	6306-6080
7	6307-6030	6307-6040	6307-6060	6307-6080
8	6308-6030	6308-6040	6308-6060	6308-6080
9	6309-6030	6309-6040	6309-6060	6309-6080
10	6310-6030	6310-6040	6310-6060	6310-6080
12	6312-6030	6312-6040	6312-6060	-
14	-	6314-6040	6314-6060	6314-6080

ORDER CODES - 120CM APPLICATION DEVICE

Stent Ø (mm)	Stent length (mm)			
	30	40	60	80
6	6306-7030	6306-7040	6306-7060	6306-7080
7	6307-7030	6307-7040	6307-7060	6307-7080
8	6308-7030	6308-7040	6308-7060	6308-7080
9	6309-7030	6309-7040	6309-7060	6309-7080
10	6310-7030	6310-7040	6310-7060	6310-7080
12	6312-7030	6312-7040	6312-7060	-
14	-	6314-7040	6314-7060	6314-7080



6F / 75 & 120cm
delivery system





Adapted to
0.035 inch
guide wire



Box / 1 unit

**FOR MORE VASCULAR INTERVENTION
HIGHLIGHTS VISIT OUR WEBSITE**

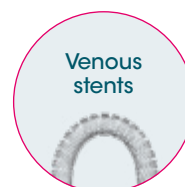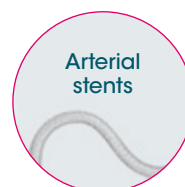

optimed Medizinische Instrumente GmbH
Ferdinand-Porsche-Straße 11 | 76275 Ettlingen Germany
Tel. +49 (0) 7243 / 76 33-0 | Fax +49 (0) 7243 / 76 33-624
info@optimed.com | www.optimed.com

CE 0124

Order Code 9997-0021EN
Rev. 08/2024 | © optimed

